YCProductSDK 修订记录

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
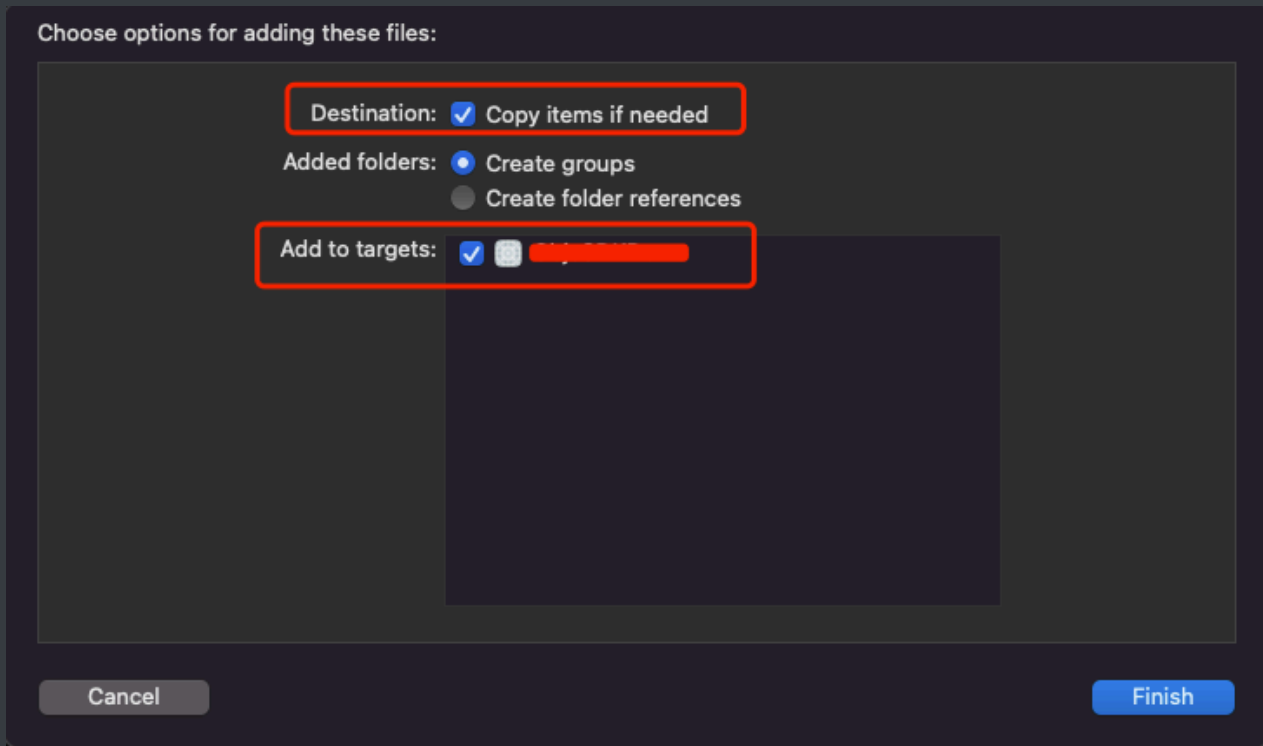
```
platform :ios, '9.0'

target 'Your app' do
  use_frameworks!

  pod 'YCPProductSDK-Swift'

end
```

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```
NSBluetoothAlwaysUsageDescription (iOS 13+)
NSBluetoothPeripheralUsageDescription
```

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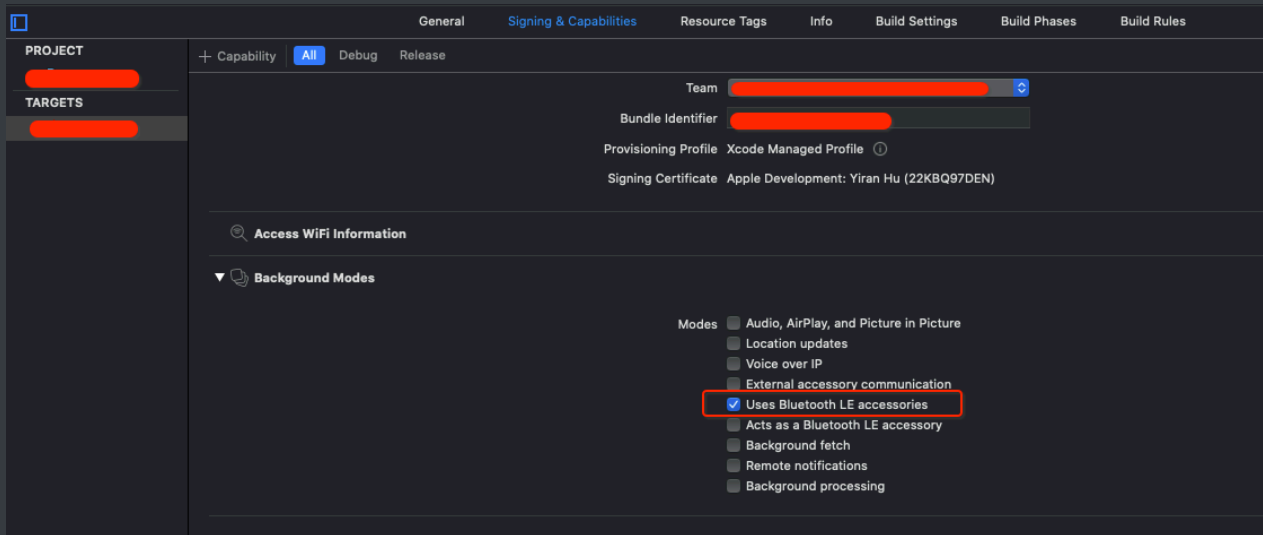```

YCProduct.scanningDevice(delayTime: 3.0) { devices, error in
    for device in devices {
        print(device.name ?? "", device.identifier.uuidString)
    }
}

```

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```
/// 连接设备
/// - Parameter peripheral: 用户选择需要连接的设备
public static func connectDevice(
    _ peripheral: CBPeripheral,
    completion: ((YCProductSDK.YCProductState, NSError?) -> ())?
)
```

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```
YCProduct.disconnectDevice(devcie) { state, _ in  
  
}
```

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```
/// 当前的连接的外设  
public var currentPeripheral: CBPeripheral? { get }  
  
/// 当前连接的所有的外设（保留参数，扩展使用）  
public var connectedPeripherals: [CBPeripheral] { get }
```

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```

        _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
        dataType: YCQueryHealthDataType,
        completion: ((_ state: YCProductState, _ response: Any?) -> (()))?
    )

```

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```

    if state == .succeed, let datas = response as? [YCHealthDataSleep] {
        for info in datas {
            print(info.startTimeStamp,
                  info.endTimeStamp,
                  info.lightSleepCount,
                  info.lightSleepMinutes,
                  info.deepSleepCount,
                  info.deepSleepMinutes,
                  info.sleepDetailDatas
                )
        }
    }
}

```

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```

YCProduct.queryHealthData(dataType: YCQueryHealthDataType.combinedData) {
state, response in

```

```

    if state == .succeed, let datas = response as?
[YCHealthDataCombinedData] {
        for info in datas {
            print(info.startTimeStamp,
                  info.bloodOxygen,
                  info.respirationRate,
                  info.temperature,
                  info.fat
                )
        }
    }
}

```

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```
YCProduct.queryDeviceBasicInfo { state, response in
    if state == YCProductState.succeed,
        let info = response as? YCDeviceBasicInfo {
            print(info.batteryPower)
        }
}
```

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```
/// 获取mac地址
public static func queryDeviceMacAddress(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((YCProductSDK.YCProductState, Any?) -> ()))?
)

/// mac地址属性
extension CBPeripheral {
    public var macAddress: String
}
```

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```
YCProduct.queryDeviceMacAddress { state, response in
    if state == YCProductState.succeed,
        let macaddress = response as? String {
            print(macaddress)
        }
}
```

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```

/// 获取型号
public static func queryDeviceModel(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((YCProductSDK.YCProductState, Any?) -> ())?
)

```

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```

YCProduct.queryDeviceModel { state, response in
    if state == YCProductState.succeed,
        let name = response as? String {
            print(name)
        }
}

```

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```

/// 获取主题
public static func queryDeviceTheme(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

/// 主题信息
@objcMembers public class YCDeviceTheme : NSObject {

    /// 总数
    public var themeCount: Int { get }

    /// 当前主题索引
    public var themeIndex: Int { get }
}

```

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```
YCProduct.queryDeviceTheme { state, response in
    if state == YCProductState.succeed,
        let info = response as? YCDeviceTheme {
            print(info.themeCount, info.themeIndex)
        }
}
```

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```
/// 获取芯片型号
public static func queryDeviceMCU(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// MCU
@objc public enum YCDeviceMCUType : UInt8 {
    case nrf52832
    case rtk8762c
    case rtk8762d
    case jl701n
}
```

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```

YCProduct.queryDeviceHistorySummary { state, response in
    if state == .succeed,
    let info = response as? YCDeviceHistorySummary {
        print(info.sleepTime, info.sportCount, info.heartRateCount)
    }
}

```

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```

YCProduct.queryDeviceCurrentExerciseInfo { state, response in
    if state == .succeed,
    let info = response as? YCDeviceCurrentExercise {
        print(info.step, info.calories, info.distance)
    }
}

```

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```

/// 获取当前心率
/// - Parameters:
///   - peripheral: 连接的设备
///   - completion: 结果

```

```

public static func queryDeviceCurrentHeartRate(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 当前心率
@objcMembers public class YCDeviceCurrentHeartRate : NSObject {

    /// 是否正在测量
    public var isMeasuring: Bool { get }

    /// 心率值
    public var heartRate: UInt8 { get }
}

```

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```

YCProduct.queryDeviceCurrentHeartRate { state, response in
    if state == YCProductState.succeed,
        let info = response as? YCDeviceCurrentHeartRate {
            print(info.isMeasuring, info.heartRate)
        }
}

```

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```

/// 获取当前血压
/// - Parameters:
///   - peripheral: 连接的设备
///   - completion: 结果
public static func queryDeviceCurrentBloodPressure(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 血压
@objcMembers public class YCDeviceCurrentBloodPressure : NSObject {

```

```

    /// 是否正在测量
    public var isMeasuring: Bool { get }

    /// 收缩压
    public var systolicBloodPressure: UInt8 { get }

    /// 舒张压
    public var diastolicBloodPressure: UInt8 { get }
}

```

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```

YCProduct.queryDeviceBloodOxygen { state, response in
    if state == .succeed,
        let info = response as? YCDeviceBloodOxygen {
            print(info.bloodOxygen)
        }
}

```

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```

/// 获取当前环境光
/// - Parameters:
///   - peripheral: 已连接设备
///   - completion: 结果
public static func queryDeviceAmbientLight(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 环境光
@objcMembers public class YCDeviceAmbientLight : NSObject {

    /// 是否正在测量
    public var isMeasuring: Bool { get }

    /// 环境光
    public var ambientLight: UInt16 { get }
}

```

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```

        public var acquisitionInterval: UInt16 { get }
    }

```

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```
YCProduct.queryDeviceWorkMode { state, response in
    if state == .succeed,
        let info = response as? YCDeviceWorkModeType {
            print(info)
        }
}
```

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```
/// 获取上传提醒配置信息
/// - Parameters:
///   - peripheral: 已连接设备
///   - completion: 结果
public static func queryDeviceUploadReminderInfo(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 上传提醒配
@objcMembers public class YCDeviceUploadReminderInfo {

    /// 是否开启
    public var isOn: Bool { get }

    /// 存储阈值 0 ~ 100
    public var threshold: UInt8 { get }
}
```

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```
YCProduct.queryDeviceRemindSettingInfo(dataType: .deviceDisconnected) {
    state, response in
        if state == .succeed,
            let state = response as? YCDeviceReminderSettingState {
                print(state == .on)
            }
}
```

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```
/// 获取屏幕显示信息
/// - Parameters:
///   - peripheral: 已连接设备
///   - completion: 结果
public static func queryDeviceScreenDisplayInfo(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 屏幕显示信息
@objc public class YCDeviceScreenDisplayInfo : NSObject {

    /// 显示等级
    public var brightnessLevel:
        YCProductSDK.YCDeviceDisplayBrightnessLevel { get }

    /// 息屏时间
    public var restScreenTime: YCProductSDK.YCDeviceBreathScreenInterval {
        get }

    /// 语言
    public var language: YCProductSDK.YCDeviceLanguageType { get }

    /// 工作模式
    public var workmode: YCProductSDK.YCDeviceWorkModeType { get }
}
```

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```
YCProduct.queryDeviceScreenDisplayInfo { state, response in
    if state == .succeed,
        let info = response as? YCDeviceScreenDisplayInfo {
            print(info.brightnessLevel, info.restScreenTime, info.language,
                info.workmode)
        }
}
```

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```
/// 时间设置
/// - Parameters:
///   - peripheral: 外设
///   - year: 年 2000+
///   - month: 月 1 ~ 12
///   - day: 日 1 ~ 31
///   - hour: 时(24) 0 ~ 23
///   - minute: 分 0 ~ 59
///   - second: 分 0 ~ 59
///   - weekDay: 星期
///   - completion: 设置是否成功
public static func setTime(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    year: UInt16,
    month: UInt8,
    day: UInt8,
    hour: UInt8,
    minute: UInt8,
    second: UInt8,
    weekDay: YCWeekDay,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

/// 设置星期
@objc public enum YCWeekDay : UInt8 {
    case monday
    case tuesday
    case wednesday
    case thursday
    case friday
    case saturday
    case sunday
}
```

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```
// 2021/12/6 14:38:59 星期一
YCProduct.setDeviceTime(
    year: 2021,
    month: 12,
    day: 6,
    hour: 14,
    minute: 38,
    second: 59,
    weekDay: .monday) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
```

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```
/// 设置步数目标
/// - Parameters:
///   - peripheral: 已连接设备
///   - step: 步数目标(步)
///   - completion: 设置结果
public static func setDeviceStepGoal(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    step: UInt32,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 设置卡路里目标
/// - Parameters:
///   - peripheral: 已连接设备
///   - calories: 卡路里目标(千卡)
///   - completion: 设置结果
public static func setDeviceCaloriesGoal(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    calories: UInt32,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
```

```

)

/// 设置距离目标
/// - Parameters:
///   - peripheral: 已连接设备
///   - calories: 距离目标(米)
///   - completion: 设置结果
public static func setDeviceDistanceGoal(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    distance: UInt32,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 设置睡眠目标
/// - Parameters:
///   - peripheral: 已连接设备
///   - hour: 时 0~ 23
///   - minute: 分 0 ~ 59
///   - completion: 设置结果
public static func setDeviceSleepGoal(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    hour: UInt8,
    minute: UInt8,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 运动时间目标
/// - Parameters:
///   - peripheral: 已连接设备
///   - hour: 运动时间 小时部分
///   - minute: 运动时间 分钟部分
///   - completion: 设置结果
public static func setDeviceSportTimeGoal(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    hour: UInt8,
    minute: UInt8,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 有效步数目标
/// - Parameters:
///   - peripheral: 已连接设备
///   - effectiveSteps: 有效步数目标 (步)
///   - completion: 设置结果

```

```
public static func setDeviceEffectiveStepsGoal(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    effectiveSteps: UInt32,
    completion: ((_ state: YCProductState, _ response: Any?) -> (()))?
)
```

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```

        completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
    )

    /// 性别
    @objc public enum YCDeviceGender: UInt8 {
        case male
        case female
    }

```

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```

YCProduct.setDeviceSedentary(startHour1: 9,
    startMinute1: 0,
    endHour1: 12,
    endMinute1: 30,
    startHour2: 13,
    startMinute2: 30,
    endHour2: 18,
    endMinute2: 00,
    interval: 15,
    repeat: [
        .monday,
        .tuesday,
        .wednesday,
        .thursday,
        .friday,
        .enable
    ]
) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}

```

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```

/// - maxHeartRate: 心率告警上限 100 ~ 240
/// - minHeartRate: 心率告警下限 30 ~ 60
/// - completion: 设置结果
public static func setDeviceHeartRateAlarm(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    isEnabled: Bool,
    maxHeartRate: UInt8,
    minHeartRate: UInt8 ,
    completion: ((_ state: YCProductState, _ response: Any?) -> (()))?
)

```

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```

highTemperatureIntegerValue: UInt8,
highTemperatureDecimalValue: UInt8,
lowTemperatureIntegerValue: Int8,
lowTemperatureDecimalValue: UInt8,
completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

```

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```
/// 恢复出厂设置
/// - Parameters:
///   - peripheral: 已连接设备
///   - completion: 结果
public static func setDeviceReset(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)
```

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```
/// 设置传感器采样率
/// - Parameters:
///   - peripheral: 已连接设备
///   - ppg: PPG采样率 HZ
///   - ecg: ECG采样率 HZ
///   - gSensor: G-Sensor采样率 HZ
///   - tempeatureSensor: 温度传感器采样率 HZ
///   - completion: 设置结果
public static func setDeviceSensorSamplingRate(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    ppg: UInt16,
    ecg: UInt16,
    gSensor: UInt16,
    tempeatureSensor: UInt16,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)
```

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```

YCProduct.setDeviceSleepReminder(hour: 22,
                                   minute: 30,
                                   repeat: [.monday, .thursday, .wednesday,
                                           .enable]) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}

```

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```
/// 查询闹钟
/// - Parameters:
///   - peripheral: 已连接设备
///   - completion: 闹钟信息
public static func queryDeviceAlarmInfo(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)
```

```

/// 闹钟信息
@objc public class YCDeviceAlarmInfo : NSObject {

    /// 设备允许的最多闹钟数量
    public var limitCount: UInt8 { get }

    /// 闹钟类型
    public var alarmType: YCProductSDK.YCDeviceAlarmType { get }

    /// 闹钟中的小时（24小时制）
    public var hour: UInt8 { get }

    /// 闹钟中的分钟
    public var minute: UInt8 { get }

    /// 重复
    public var `repeat`: Set<YCProductSDK.YCDeviceWeekRepeat> { get }

    /// 贪睡时间 min
    public var snoozeTime: UInt8 { get }
}

/// 闹钟类型
@objc public enum YCDeviceAlarmType : UInt8 {
    case wakeUp           // 起床
    case sleep             // 睡眠
    case exercise          // 锻炼
    case medicine          // 吃药
    case appointment      // 约会
    case party             // 聚会
    case meeting           // 会议
    case custom            // 自定义
}

```

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```

YCPProduct.addDeviceAlarm(alarmType: .wakeUp,
                           hour: 6,
                           minute: 30,
                           repeat: [.enable, .sunday, .saturday],
                           snoozeTime: 0) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail \(String(describing: response) ?? "")")
    }
}

```

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```

/// 修改闹钟
/// - Parameters:
///   - peripheral: 已连接设备
///   - oldHour: 闹钟原来的小时
///   - oldMinute: 闹钟原来的分钟
///   - hour: 闹钟新的小时
///   - minute: 闹钟新的分钟
///   - alarmType: 闹钟类型
///   - repeat: 重复时间
///   - snoozeTime: 贪睡时长
///   - completion: 结果
public static func modifyDeviceAlarm(
    _ peripheral: CBPeripheral? = YCPProduct.shared.currentPeripheral,
    oldHour: UInt8,
    oldMinute: UInt8,
    hour: UInt8,
    minute: UInt8,
    alarmType: YCDeviceAlarmType,
    repeat: Set<YCDeviceWeekRepeat>,
    snoozeTime: UInt8,
    completion: ((_ state: YCPProductState, _ response: Any?) -> ()))?
)

```

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```
YCProduct.deleteDeviceAlarm(hour: 6, minute: 30) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
```

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```

///    - minimum: 报警值下限
///    - completion: 结果
public static func setDeviceRespirationRateAlarm(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    isEnabled: Bool,
    maximum: UInt8,
    minimum: UInt8,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

```

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```

YCProduct.queryDeviceRealTimeTemperature { state, response in
    if state == .succeed,
        let temperature = response as? Double {
        print("success \((temperature)")
    } else {
        print("fail")
    }
}

```

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```

YCProduct.deviceSystemOperator(mode: .shutDown) { state, _ in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
}

```

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```

/// 接收实时数据的通知
public static let receivedRealTimeNotification: Notification.Name

/// 接收到实时数据类型
@objc public enum YCReceivedRealTimeDataType: UInt8 {
    case step
    case heartRate
    case bloodOxygen
    case bloodPressure
    case ppg
    case ecg
    case comprehensiveData
    case realTimeMonitoringMode
}

/// 接收到的实时数据的响应
@objcMembers public class YCReceivedDeviceInfo: NSObject {

    /// 收到的响应设备
    public var device: CBPeripheral? { get }

    /// 结果
    public var data: Any? { get }
}

/// 实时步数
@objcMembers public class YCReceivedRealTimeStepInfo: NSObject {

    /// 步数（单位:步）

```

```

    /// 血氧值
    public var bloodOxygen: Int { get }

    /// 呼吸率值
    public var respirationRate: Int { get }

    /// 温度
    public var temperature: Double { get }

    /// 实时步数
    public var realStep: Int { get }

    /// 实时距离 (单位:米)
    public var realDistance: UInt16 { get }

    /// 实时卡路里 (单位:千卡)
    public var realCalories: UInt16 { get }

    /// 模式步数
    public var modeStep: Int { get }

    /// 模式距离 (单位:米)
    public var modeDistance: UInt16 { get }

    /// 模式卡路里 (单位:千卡)
    public var modeCalories: UInt16 { get }

    /// 脉搏波信号的峰峰值间期 (单位: 微秒)
    public var ppi: Int { get }
}

```

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```
NotificationCenter.default.addObserver(
```

```

        self,
        selector: #selector(receiveRealTimeData(_:)),
        name: YCProduct.receivedRealTimeNotification,
        object: nil
    )

@objc private func receiveRealTimeData(_ notification: Notification) {

    guard let info = notification.userInfo else {
        return
    }

    if let response = info[YCReceivedRealTimeDataType.step.string] as?
YCReceivedDeviceInfo,
        let device = response.device,
        let sportInfo = response.data as? YCReceivedRealTimeStepInfo {

        print(device.name ?? "",
            sportInfo.step,
            sportInfo.calories,
            sportInfo.distance
        )
    }

    else if let response =
info[YCReceivedRealTimeDataType.heartRate.string] as?
YCReceivedDeviceInfo,
        let device = response.device,
        let heartRate = response.data as? UInt8 {

        print(device.name ?? "",
            heartRate)
    }

    else if let response =
info[YCReceivedRealTimeDataType.bloodOxygen.string] as?
YCReceivedDeviceInfo,
        let device = response.device,
        let bloodOxygen = response.data as? UInt8 {

        print(device.name ?? "",
            bloodOxygen)
    }
}

```

```

        else if let response =
            info[YCReceivedRealTimeDataType.bloodPressure.string] as?
            YCReceivedDeviceInfo,
            let device = response.device,
            let bloodPressureInfo = response.data as?
            YCReceivedRealTimeBloodPressureInfo {

                print(device.name ?? "",
                    bloodPressureInfo.systolicBloodPressure,
                    bloodPressureInfo.diastolicBloodPressure)
            }
        }
    }
}

```

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```

/// 波形上传控制
/// - Parameters:
///   - peripheral: 已连接设备
///   - state: 是否开关
///   - dataType: 波形类型
///   - completion: 结果
public static func waveDataUpload(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    state: YCWaveUploadState,
    dataType: YCWaveDataType,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

/// 波形类型选择
@objc public enum YCWaveDataType: UInt8 {
    case ppg
    case ecg
    case multiAxisSensor
    case ambientLight
}

@objc public enum YCWaveUploadState: UInt8 {
    case off = 0 // 停止传输
    case uploadWithoutSerialnumber // 无序号传输
    case uploadSerialnumber // 带8位充号传输
}

```

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```
NotificationCenter.default.addObserver(
    self,
    selector: #selector(receiveRealTimeData(_:)),
    name: YCProduct.receivedRealTimeNotification,
    object: nil
)
```

```
@objc private func receiveRealTimeData(_ notification: Notification) {

    guard let info = notification.userInfo else {
        return
    }

    if let response = info[YCReceivedRealTimeDataType.ppg.string] as?
YCReceivedDeviceInfo,
        let device = response.device,
```

```

        let ppgData = response.data as? [Int32] {
            print(device.name ?? "", ppgData)

        } else if let response = info[YCReceivedRealTimeDataType.ecg.string]
as? YCReceivedDeviceReportInfo,
            let device = response.device,
            let ecgData = response.data as? [Int32] {
                print(device.name ?? "", ecgData)
            }
    }
}

```

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```
// 左手
YCProduct.setDeviceWearingPosition(wearingPosition: .left) { state,
response in
    if state == .succeed {

    }
}
```

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```

/// 全局对象
public static let shared: YCProductSDK.YCECGManager

/// 算法计算过程回调
public func setupManagerInfo(
    rr:((_ rr: Float, _ heartRate: Int) -> ()))?,
    hrv: ((_ hrv: Int) -> ())?
)

```

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```
/// 接收实时数据的通知
public static let receivedRealTimeNotification: Notification.Name
```

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```
NotificationCenter.default.addObserver(
    self,
    selector: #selector(receiveRealTimeData(_:)),
    name: YCProduct.receivedRealTimeNotification,
    object: nil
)

@objc private func receiveRealTimeData(_ notification: Notification) {
    guard let info = notification.userInfo else {
        return
    }
}
```

```

    }
    if let response = info[YCReceivedRealTimeDataType.step.string] as?
YCReceivedDeviceInfo,
        let device = response.device,
        let sportInfo = response.data as? YCReceivedRealTimeStepInfo {
            print(device.name ?? "",
                sportInfo.step,
                sportInfo.calories,
                sportInfo.distance
            )
        }
    else if let response =
info[YCReceivedRealTimeDataType.heartRate.string] as?
YCReceivedDeviceInfo,
        let device = response.device,
        let heartRate = response.data as? UInt8 {
            print(device.name ?? "",
                heartRate)
        }
    }
}

```

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```

/// 接收实时数据的通知
public static let receivedRealTimeNotification: Notification.Name

/// 运动状态变化通知
public static let deviceControlNotification: Notification.Name

/// 运动状态信息
@objcMembers public class YCDeviceControlSportModeControlInfo: NSObject {
    public var state: YCDeviceSportState {get}
    public var sportType: YCDeviceSportType {get}
}

```

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```

NotificationCenter.default.addObserver(
    self,

```

```

        selector: #selector(deviceDataStateChanged(_:)),
        name: YCProduct.deviceControlNotification,
        object: nil
    )

```

```

NotificationCenter.default.addObserver(
    self,
    selector: #selector(receiveRealTimeData(_:)),
    name: YCProduct.receivedRealTimeNotification,
    object: nil
)

```

```

@objc private func receiveRealTimeData(_ notification: Notification) {

```

```

    guard let info = notification.userInfo else {
        return
    }

```

```

    if let response =
info[YCReceivedRealTimeDataType.realTimeMonitoringMode.string] as?
YCReceivedDeviceReportInfo,
        let device = response.device,
        let data = response.data as? YCReceivedMonitoringModeInfo {
        print(device.name ?? "",
            data.startTimeStamp,
            data.modeStep,
            data.modeCalories,
            data.modeCalories
        )
    }
}

```

```

@objc private func deviceDataStateChanged(_ ntf: Notification) {

```

```

    guard let info = ntf.userInfo else {
        return
    }

```

```

    if let response = info[YCDeviceControlType.sportModeControl.string] as?
YCReceivedDeviceReportInfo,
        let device = response.device,
        let data = response.data as? YCDeviceControlSportModeControlInfo {
        print(device.name ?? "",
            data.state,

```

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```

/// 手机控制进入与退出拍照模式
/// - Parameters:
///   - peripheral: 已连接设备
///   - isEnabled: 启动或关闭拍照模式
///   - completion: 结果
public static func takephotoByPhone(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    isEnabled: Bool,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

```

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```
/// 发送健康参数
/// - Parameters:
///   - peripheral: 已连接设备
///   - warningState: 预警状态
///   - healthState: 健康状态
///   - healthIndex: 健康指数 0 ~ 120
///   - othersWarningState: 亲友预警是否生效
///   - completion: 结果
public static func sendHealthParameters(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    warningState: YCHealthParametersState,
    healthState: YCHealthState,
    healthIndex: UInt8,
    othersWarningState: YCHealthParametersState,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 预警状态
@objc public enum YCHealthParametersState: UInt8 {
    case off          // 无预警
    case effect        // 预警生效中
    case invalid       // 无效
}

/// 健康状态
@objc public enum YCHealthState: UInt8 {
    case unknow        // 未知
    case excellent      // 优秀
    case good           // 良好
    case general        // 一般
    case poor           // 较差
    case sick           // 生病
    case invalid        // 无效
}
```

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```

YCProduct.sendHealthParameters(warningState: .off, healthState: .good,
healthIndex: 100, othersWarningState: .off) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
}

```

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```

/// 显示亲友消息
/// - Parameters:
///   - peripheral: 已连接设备
///   - index: 表情序号 0 ~ 4
///   - hour: 发送小时 0 ~ 23
///   - minute: 发送分钟 0 ~ 59
///   - name: 亲友名称
///   - completion: 结果
public static func deviceShowFriendMessage(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,

    index: UInt8,
    hour: UInt8,
    minute: UInt8,
    name: String,
    completion: ((_ state: YCProductState, _ response: Any?) -> (()))?
)

```

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```

YCProduct.deviceShowFriendMessage(index: 1,
                                   hour: 10,
                                   minute: 23,
                                   name: "俺是大宝二") { state, response in

    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}

```

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```

YCProduct.deviceHealthValueWriteBack(healthValue: 50, statusDescription:
"非常好") { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
}

```

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```

/// 睡眠数据回写
/// - Parameters:
///   - peripheral: 已连接设备
///   - deepSleepHour: 深睡时长 小时 0 ~ 23
///   - deepSleepMinute: 深睡时长 分钟 0 ~ 23
///   - lightSleepHour: 浅睡时长 小时 0 ~ 23
///   - lightSleepMinute: 浅睡时长 分钟 0 ~ 59
///   - totalSleepHour: 浅睡时长 小时 0 ~ 23
///   - totalSleepMinute: 总共睡眠时长 分钟 0 ~ 59
///   - completion: 结果
public static func deviceSleepDataWriteBack(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,

    deepSleepHour: UInt8,
    deepSleepMinute: UInt8,
    lightSleepHour: UInt8,
    lightSleepMinute: UInt8,
    totalSleepHour: UInt8,
    totalSleepMinute: UInt8,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

```

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```

YCProduct.deviceSleepDataWriteBack(
    deepSleepHour: 2,
    deepSleepMinute: 30,
    lightSleepHour: 4,
    lightSleepMinute: 0,
    totalSleepHour: 8,
    totalSleepMinute: 0) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}

```

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```

/// 个人数据回写
/// - Parameters:
///   - peripheral: 已连接设备
///   - infoType: 用户信息类型
///   - information: 描述信息
///   - completion: 结果
public static func devicePersonalInfoWriteBack(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,

    infoType: YCPersonalInfoType,

    information: String,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

/// 个人信息类型
@objc public enum YCPersonalInfoType: UInt8 {
    case insurance          // 保险
    case vip                // vip
}

```

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```

/// 个人数据回写
/// - Parameters:
///   - peripheral: 已连接设备
///   - infoType: 用户信息
///   - information: 描述信息
///   - completion: 结果
public static func devicePersonalInfoWriteBack(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,

    infoType: YCPersonalInfoType,
    information: String,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

```

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```

/// 升级提醒
/// - Parameters:
///   - peripheral: 已连接设备
///   - isEnabled: 是否开启提醒
///   - percentage: 当前进度 0 ~ 100
///   - completion: 结果
public static func deviceUpgradeReminderWriteBack(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,

    isEnabled: Bool,
    percentage: UInt8,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

```

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```

YCProduct.deviceUpgradeReminderWriteBack(isEnable: true,
                                         percentage: 60) { state, response
in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}

```

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```

/// 运动数据回写
/// - Parameters:
///   - peripheral: 已连接设备
///   - step: 运动步数
///   - state: 运动状态
///   - completion: 结果
public static func deviceSportDataWriteBack(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    step: UInt32,
    state: YCDeviceExerciseHeartRateType,

    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

/// 运动状态类型
@objc public enum YCDeviceExerciseHeartRateType: UInt8 {
    case retreat // 休养静
歇
    case casualwarmup // 休闲热
身
    case cardiorespiratoryStrengthening // 心肺强化
    case reduceFatShape // 减脂塑形

    case sportsLimit // 运
动极限
    case emptyState // 空状
态
}

```

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```
YCProduct.deviceSportDataWriteBack(step: 10000, state: .reduceFatShape) {
state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
```

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```
/// 计算心率发送到设备
/// - Parameters:
///   - peripheral: 已连接设备
///   - heartRate: 计算心率
///   - completion: 结果
public static func sendCaclulateHeartRate(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    heartRate: UInt8,
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)
```

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```
YCProduct.sendCaclulateHeartRate(heartRate: 78) { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
```

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```
/// 健康数据回写
/// - Parameters:
///   - peripheral: 已连接设备
///   - dataType: 测量数据类型
///   - values: 测量值集合
///   - completion: 结果
public static func deviceMeasurementDataWriteBack(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    dataType: YCMeasurementDataType,
    values: [UInt8],
    completion: ((_ state: YCProductState, _ response: Any?) -> ())?
)

/// 健康数据测量类型
@objc public enum YCMeasurementDataType: UInt8 {

    case heartRate          = 0    // 心率
    case bloodPressure      // 血压
    case bloodOxygen        // 血氧
    case respirationRate    // 呼吸率
    case hrv                 // hrv
    case bloodGlucose        // 血糖
    case temperature        // 温度
}
```

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```
/// 发送手机型号
/// - Parameters:
///   - peripheral: 已连接设备
///   - mode: 手机型号 如iPhone 13
///   - completion: 结果
public static func sendPhoneModeInfo(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    mode: String,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)
```

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```

YCProduct.sendPhoneModeInfo(mode: "iPhone13 Pro Max") { state, response in
    if state == .succeed {
        print("success")
    } else {
        print("fail")
    }
}
}

```

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```

/// 预警信息
/// - Parameters:
///   - peripheral: 已连接设备
///   - infoType: 信息类型
///   - message: 信息内容
///   - completion: 结果
public static func sendWarningInformation(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    infoType: YCWarningInformationType,
    message: String?,
    completion: ((_ state: YCProductState, _ response: Any?) -> ()))?
)

/// 预警信息类型
@objc public enum YCWarningInformationType: UInt8 {
    case warnSelf          // 预警自己
    case warnOthers        // 预警他人
    case highRisk           // 运动高风险
    case nonHighRisk       // 非高风险
}

```

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```

/// 开启通讯录同步
/// - Parameters:
///   - peripheral: 已连接设备
///   - completion: 结果
public static func startSendAddressBook(
    _ peripheral: CBPeripheral? = YCPProduct.shared.currentPeripheral,
    completion: ((_ state: YCPProductState, _ response: Any?) -> ())?
)

```

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```

        if let response = info[YCDeviceControlType.findPhone.toString] as?
YCReceivedDeviceInfo,
            let device = response.device,
            let state = response.data as? YCDeviceControlState {
                print(device.name ?? "",
                    state == .stop
                )
            }

// sos
        if let response = info[YCDeviceControlType.sos.toString] as?
YCReceivedDeviceInfo,
            let device = response.device {
                print(device.name ?? "",
                    "sos"
                )
            }

// 是否允许连接
        if let response = info[YCDeviceControlType.allowConnection.toString] as?
YCReceivedDeviceInfo,
            let device = response.device,
            let state = response.data as? YCDeviceControlAllowConnectionState {
                print(device.name ?? "",
                    state == .agree
                )
            }

// 恢复出厂设置 reset
        if let response = info[YCDeviceControlType.reset.toString] as?
YCReceivedDeviceInfo,
            let device = response.device {
                print(device.name ?? "",
                    "reset"
                )
            }

// 预警值
        if let response = info[YCDeviceControlType.reportWarningValue.toString]
as? YCReceivedDeviceInfo,
            let device = response.device,
            let value = response.data as? YCDeviceControlReportWarningValueInfo {
                print(device.name ?? "",
                    value
                )
            }

```

```

        )
    }

    if let response =
info[YCDeviceControlType.invasiveMeasurementState.toString] as?
YCReceivedDeviceInfo,
        let device = response.device,
        let value = response.data as?
YCDeviceControlInvasiveMeasurementStateInfo {

        print(value.toString)
    }
}

```

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```

@objc public enum YCDeviceControlAllowConnectionState: UInt8 {
    case agree
    case refuse
}

```

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```
@objcMembers public class YCDeviceControlReportWarningValueInfo: NSObject {

    /// 预警类型
    public var dataType: YCAppControlMeasureHealthDataType {get}

    /// 预警值
    public var values: [Int] {get}
}

// 如果测量类型是血压，values中的第一个元素是收缩压，第二个是舒张压。
```

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```

/// 删除表盘
/// - Parameters:
///   - peripheral: 已连接设备
///   - dialID: 表盘ID
///   - completion: 删除结果
public static func deleteWatchFace(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    dialID: UInt32,
    completion: ((YCProductSDK.YCProductState, Any?) -> ())?
)

```

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```

let dialID: UInt32 = 2147483539
YCProduct.deleteWatchFace(dialID: dialID) { state, _ in
    if state == .succeed {
        print("delete success")
    }
}

```

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```
/// 切换表盘
/// - Parameters:
///   - peripheral: 已连接设备
///   - dialID: 表盘ID
///   - completion: 结果
public static func changeWatchFace(
    _ peripheral: CBPeripheral? = YCProduct.shared.currentPeripheral,
    dialID: UInt32,
    completion: ((YCProductSDK.YCProductState, Any?) -> ())?
)
```

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```
/// 设备控制通知
public static let deviceControlNotification: Notification.Name
```

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```

let dialID: UInt32 = 2147483539

```

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```

YCProduct.deleteWatchFace(dialID: dialID) { state, response in

```

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```

    YCProduct.downloadWatchFace(
        isEnabled: true,
        data: customDialData,
        dialID: dialID,
        blockCount: 0,
        dialVersion: 1) { state, response in
        if state == .succeed,
            let info = response as? YCDownloadProgressInfo {
            print(info.downloaded, info.progress)
        } else {

        }
    }
}

```

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```
/// 生成自定义表盘数据
/// - Parameters:
///   - dialData: 原始表盘数据
///   - backgroundImage: 背景图片
///   - thumbnail: 缩略图片
///   - timePosition: 时间显示位置坐标
///   - redColor: 0 ~ 255
///   - greenColor: 0 ~ 255
///   - blueColor: 0 ~ 255
///   - isFlipColor: 是否要翻转颜色
/// - Returns: 表盘数据
public static func generateCustomDialData(
    _ dialData: Data,
    backgroundImage: UIImage?,
    thumbnail: UIImage?,
    timePosition: CGPoint,
    redColor: UInt8,
    greenColor: UInt8,
    blueColor: UInt8,
    isFlipColor: Bool
) -> Data

/// 查询表盘文件中的BMP信息
public static func queryDeviceBmpInfo(_ dialData: Data) ->
YCProductSDK.YCWatchFaceDataBmpInfo

/// 表盘中的图片信息
@objcMembers public class YCWatchFaceDataBmpInfo : NSObject {

    /// 背景图片的宽度
    public var width: Int { get }

    /// 背景图片的高度
    public var height: Int { get }

    /// 背景图片的大小(字节)
    public var size: Int { get }

    /// 背景图片的半径
    public var radius: Int { get }
```

```

    /// 缩略图片的宽度
    public var thumbnailWidth: Int { get }

    /// 缩略图片的高度
    public var thumbnailHeight: Int { get }

    /// 缩略图片的大小(字节)
    public var thumbnailSize: Int { get }

    /// 缩略图片的半径
    public var thumbnailRadius: Int { get }
}

```

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```

    public var version: String = ""

    init(name: String, dialID: String, version: String)
}

```

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```
import iOSSDFULibrary

class YCFirmwareUpgradeViewController: UIViewController {

    /// NRF升级控制
```

```

private var dfuController: DFUServiceController?

override func viewDidDisappear(_ animated: Bool) {
    super.viewDidDisappear(animated)

    _ = dfuController?.abort()
}

}

// MARK: - 固件升级
extension YCFirmwareUpgradeViewController: DFUServiceDelegate,
DFUProgressDelegate {

    /// NRF固件升级（自定义实现）
    /// - Parameters:
    ///   - filePath: 固件文件地址
    ///   - device: 需要升级的设备
    private func startNRFFirmwareUpgrade(_ filePath: URL,
                                         device: CBPeripheral) {

        guard let dfuFirmware = DFUFirmware(urlToZipFile: filePath) else {
            printLog("固件不存在")
            return
        }

        let initiator =
            DFUServiceInitiator(queue: DispatchQueue.main,
                               delegateQueue: DispatchQueue.main,
                               progressQueue: DispatchQueue.main,
                               loggerQueue: DispatchQueue.main
                                )
        initiator.delegate = self
        initiator.progressDelegate = self
        initiator.forceDfu = false
        initiator.alternativeAdvertisingNameEnabled = false
        initiator.enableUnsafeExperimentalButtonlessServiceInSecureDfu =
true

        _ = initiator.with(firmware: dfuFirmware)
        dfuController = initiator.start(target: device)
    }

    /// 状态变化

```


```
func dfuStateDidChange(to state: DFUState) {
    switch state {

        case .disconnecting:
            break

        case .connecting:
            break

        case .starting:
            break

        case .enablingDfuMode: // 进入升级状态
            break

        case .completed: // 固件升级结束
            break

        case .aborted:
            break

        default:
            break
    }

    /// 升级出错
    func dfuError(_ error: DFUError, didOccurWithMessage message: String) {

    }

    /// 升级进度
    func dfuProgressDidChange(for part: Int, outOf totalParts: Int, to
    progress: Int, currentSpeedBytesPerSecond: Double, avgSpeedBytesPerSecond:
    Double) {

    }
}
```

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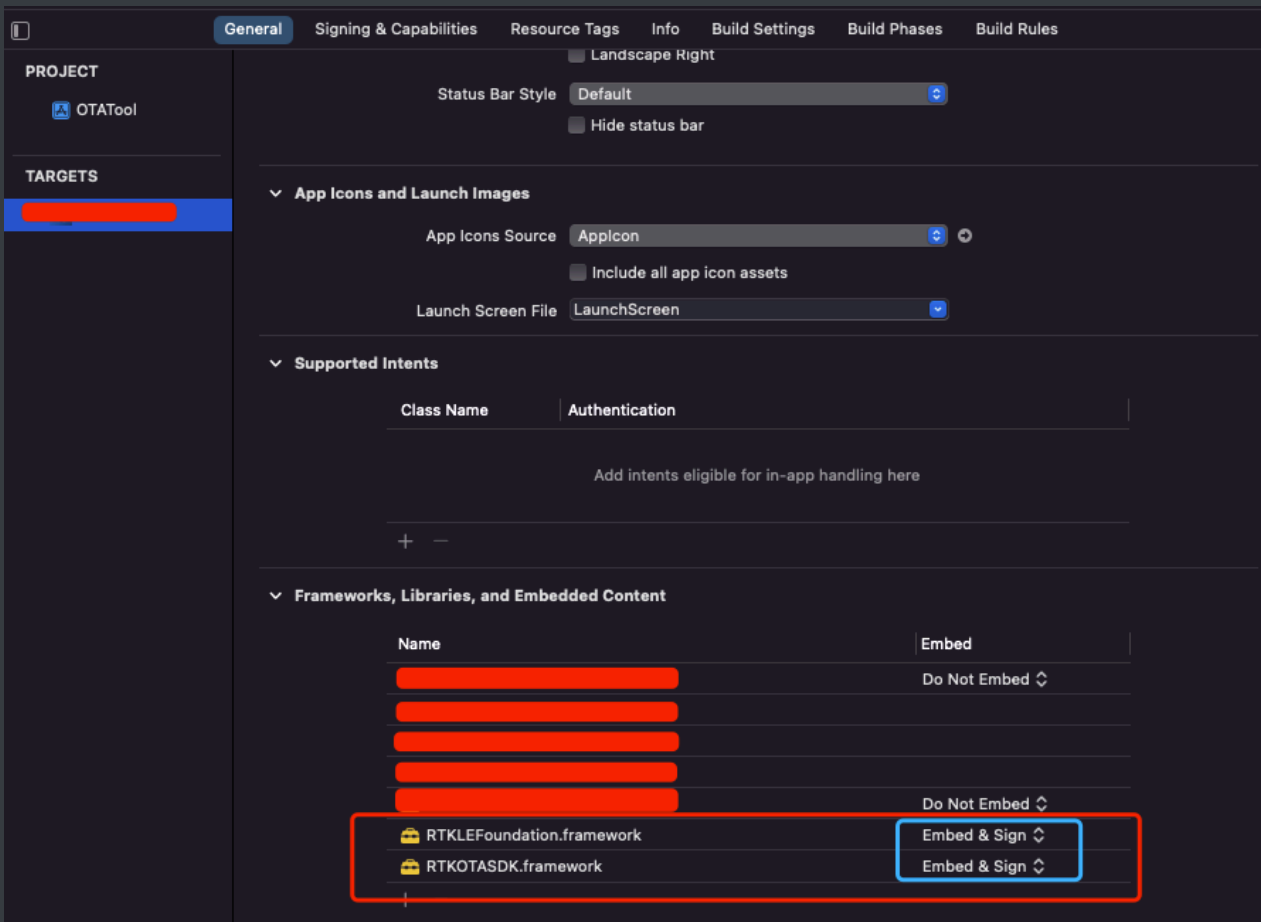```

        completion: { state, progress in

            self?.setupOTAState(
                filePath,
                state: state,
                progress: progress,
                didSend: 0
            )
        }
    )
}

// }
}

/// OTA升级
/// - Parameters:
///   - filePath: <#filePath description#>
///   - completion: <#completion description#>
private func jlDeviceOTA(_ filePath: String) {

    YCProduct.jlDeviceUpgradeFirmware(
        filePath: filePath,
        completion: { [weak self] state, progress, didSend in

            self?.setupOTAState(filePath,
                                state: state,
                                progress: progress,
                                didSend: didSend
            )
        }
    )
}

/// 升级状态
private func setupOTAState(_ filePath: String, state: JLDeviceUpdateState,
progress: Float, didSend: Float) {

    switch state {

    case .start:
        // 开始
        break
    }
}

```

```

        case .resourceUpdating:
            // 正在更新资源
            break

        case .uiUpdating:
            // 正在更新UI
            break

        case .upgrading:
            // 正在升级
            break

        case .updateUIFinished:
            // 升级UI完成，重新扫描设备，进入固件升级
            self.reconnectJLDeviceStart()

        case .failed:
            // 升级失败
            break

        case .success:
            // 升级成功
            break

        case .updateResourceFinished:
            // 更新资源结束，进入UI升级
            self.jlDeviceOTA(filePath)
            break
    }
}

// 回连设备再次升级
extension YCFirmwareUpgradeViewController {

    /// OTA回连设备
    private func reconnectJLDeviceStart() {
        // NotificationCenter.default.addObserver(
        //     self,
        //     selector: #selector(jlDevicePairedActionForOTA(_:)),
        //     name: YCProduct.jlDevicePairedNotification,
        //     object: nil
        // )
    }
}

```

```

        setupTimer()
    }

    /// 设置定时器
    private func setupTimer() {

        let timer =
            Timer(timeInterval: 8.0,
                  target: self,
                  selector: #selector(timeUp),
                  userInfo: nil,
                  repeats: true
            )

        RunLoop.current.add(timer, forMode: .common)
        stateTimer = timer
        timeCount = 0

        timer.fire()
    }

    /// 销毁定时器
    private func destoryJLTimer() {
        YCProduct.stopSearchDevice()
        stateTimer?.invalidate()
        stateTimer = nil
        timeCount = 0
    }

    /// 时间到了
    @objc private func timeUp() {
        timeCount += 8
        if timeCount >= 40 {
            // 超时，升级失败
            destoryJLTimer()
            return
        }

        // 搜索设备
        YCProduct.scanningDevice { devices, error in
            for device in devices where deviceArray.contains(device) ==
false {
                deviceArray.append(device)
            }
        }
    }

```

```

    }

    // 搜索结束 准备连接
    perform(
        #selector(connectDevice),
        with: nil,
        afterDelay: 3.5
    )
}

/// 连接设备
@objc private func connectDevice() {
    guard deviceArray.isEmpty == false,
        let reconnectDevice = reconnectDevice else {
        return
    }

    // 搜索设备
    for device in deviceArray {
        if device.macAddress.uppercased() ==
reconnectDevice.macAddress.uppercased() {

            if device.state == .connected {
                destoryJLTimer()
                break
            }

            YCProduct.connectDevice(device) { [weak self] state, error
in
                if state == .connected {
                    self?.destoryJLTimer()
                    self?.startJLFirmwareUpgrade(
                        upgradeFilePath,
                        device: device
                    )
                }
            }
            break
        }
    }
}

/*
/// OTA重新连接升级

```

```
/// - Parameter notification: <#notification description#>
@objc func jlDevicePairedActionForOTA(_ notification: Notification) {

    if let device = reconnectDevice ??
YCPProduct.shared.currentPeripheral {
        JL_Tools.delay(2.0) { [weak self] in          // 延时2秒
            NotificationCenter.default.removeObserver(self)
            self?.destroyJLTimer()
            self?.startJLFirmwareUpgrade(
                upgradeFilePath,
                device: device
            )
        }
    }
}
*/
}
```